



Button Pupilometer Model #CX-8

- 21st century ergonomic design
- Compact and lightweight
- View both monocular and binocular readings simultaneously
- Take distance PD and pupilometer will calculate near
- Memory button holds last reading until cleared
- Shuts off automatically
- Uses 4 AA batteries
- Measures PD's from 45-82mm




Specifications:

Dimensions: 8-1/2" L x 6-1/3" W x 2-1/3"H

Weight: 1.2 lbs.

- 1) **Reset Key:** press this key to turn on pupilometer.
- 2) **Right Pupil Adjustment Key:** adjusts the right pupillary distance point in .5mm increments.
- 3) **Left Pupil Adjustment Key:** adjusts the left pupillary distance point in .5mm increments.
- 4) **Display Window**
- 5) **Memory Key:** press once to store current reading and press again to clear.
- 6) **TDIS (Target Distance Selection) Key:** 30cm to ∞
- 7) **Binocular/Monocular Setting Key**
- 8) **Nose Rest**
- 9) **Eye Piece**
- 10) **Patient View Windows**

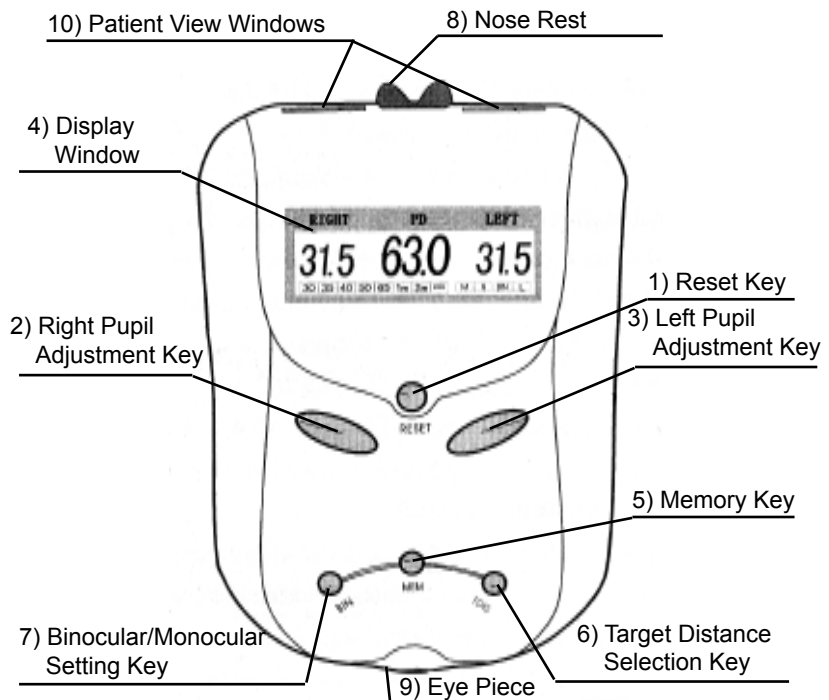

To Order Call 800-237-8154 or Fax 800-330-3800



PRACTICAL SYSTEMS INC.

11617 Prospect Road • Odessa, FL 33556
(727) 376-7900 • Fax (727) 372-1900
(800) 237-8154 • Fax (800) 330-3800
www.LookToPSI.com